





Pathology and Laboratory Medicine

Orderable - GENX

Turn Around Time: 2 weeks

Alternate Name(s):

VEGF

Specimen:

Laboratory:

Core Lab

4 mL K_2 or K_3 EDTA Lavender top Vacutainer tube

== ×-

Requisition:GENERAL LABORATORY
REQUISITION

Collection Information:

Minimum volume of plasma required is 300 μL.



Method of Analysis: Enzyme Immunoassay (EIA) **Note:** This test is available for Dr. Martha Louzada, Dr. Chai Phua, Dr. Selay Lam, Dr. Jenny Ho and Dr. Cheryl Foster. Other physicians ordering this test will require approval by the on-call biochemist prior to ordering the test.



Test Schedule: Referred out Monday-Thursday

Reference Ranges:

31 - 86 pg/mL



Referred Out Location: In-Common Laboratories

Interpretive Comments:


This test should not be ordered more frequently than 4 times per year for patients with POEMS syndrome.

This assay is intended for research use only. It may have value in the evaluation of patients with various types of cancer who are being considered for anti-angiogenesis therapy.

Specimen Stability:

Ambient: 4 hours Refrigerated: 48 hours

Frozen: 1 year

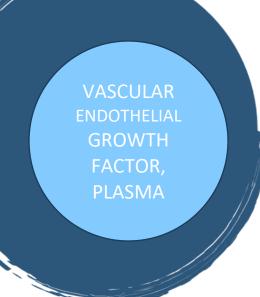





Pathology and Laboratory Medicine

Storage and Shipment:

Store and ship plasma frozen.